

# Advazorb®

## Hydrophilic foam dressing with film backing

### Product description

Advazorb® and Advazorb Lite® are soft and conformable, low-adherent, hydrophilic, polyurethane foam dressings with a breathable film backing.


A non-adherent foam dressing with excellent fluid handling properties. Ideal for use under compression bandaging with its fluid retention, low profile and low friction film backing.

The Advazorb® range features a waterproof film backing that prevents strike-through. The backing has been designed to reduce friction on clothes and bandages which extends dressing wear time by reducing 'rucking' and preventing dressings falling off.

With an MVTR range of 2600-2700 g/m<sup>2</sup>/24hours the Advazorb® foam range provides the ideal conditions for wound healing.

### Use

Apply directly to the wound surface pink film side up and secure in place with tape, appropriate bandage or secondary dressing.

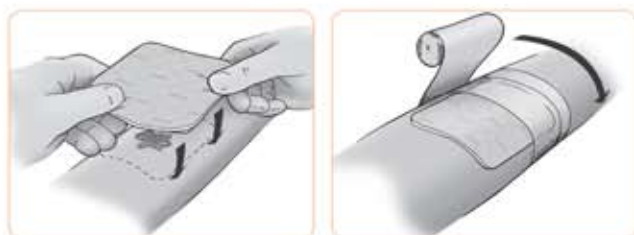

### Contra-indications

Arterial bleeds, heavily bleeding wounds and vascular fungating tumours.

### Storage

Store at room temperature, out of direct sunlight and keep dry.

### Advazorb® ordering information

Size	Stock code	Pack size	NHS code	PIP code
5cm x 5cm	CR4221	10	ELA624	374-4638
7.5cm x 7.5cm	CR4165	10	ELAS26	367-0999
10cm x 10cm	CR4166	10	ELA453	367-1005
12.5cm x 12.5cm	CR4167	10	ELAS27	367-1013
15cm x 15cm	CR4168	10	ELA454	367-1021
10cm x 10cm	CR4169	10	ELAS28	367-1039
20cm x 20cm	CR4170	10	ELA455	367-1047



#### TOP TIP

Advazorb® is ideal for use under compression bandaging with its fluid retention and low friction backing.



### Features

- ✓ Sterile Hydrophilic PU Foam
- ✓ Low friction thin film backing
- ✓ Bacteria proof backing
- ✓ Available in regular and 'Lite' versions
- ✓ Soft and conformable
- ✓ High fluid retention

### Indications

Suitable for acute and chronic exuding wounds including:

- ✓ Cuts and abrasions
- ✓ Superficial burns
- ✓ Surgical wounds
- ✓ Leg ulcers
- ✓ Pressure ulcers
- ✓ Diabetic ulcers

### Care and use symbols



**Advancis**  
Medical  
[www.advancis.co.uk](http://www.advancis.co.uk)

Another quality product selected and stocked by:

WM Bamford & Co Ltd p 0800-BAMFORD f 04-569-6489  
e enquiries@bamford.co.nz www.bamford.co.nz

When performance counts **BAMFORD** delivers

