



SILICONE FOLEY BALLOON CATHETER

The Right Choice for Quality Catheters



UROLOGICAL
SERIES



UROLOGICAL SERIES

SILICONE FOLEY BALLOON CATHETER

Fortune® Silicone 2 way Foley balloon Catheter Features:

- 100% medical grade silicone for superior biocompatibility.
- Transparent medical grade silicone tube allows easy visual inspection and fluid observation.
- X-ray opaque line allows for confirmation of intubated tube using X-ray.
- Soft and uniformly inflated balloon makes the tube sit well against the bladder.
- Smooth round shaft can minimize trauma during insertion and withdrawal.
- 6FR, 8FR and 10FR Pediatric Foley Catheter has disposable stylet for easy insertion.
- Includes an individual sterilized packed Catheter Spigot.
- Sterilized double packaging.

- A

Pediatric
- B


Standard, 5-10 cc/ml
- C

Standard, 20-30 cc/ml
- D

Female
- E

Catheter Spigot
- F


Tiemann Tip

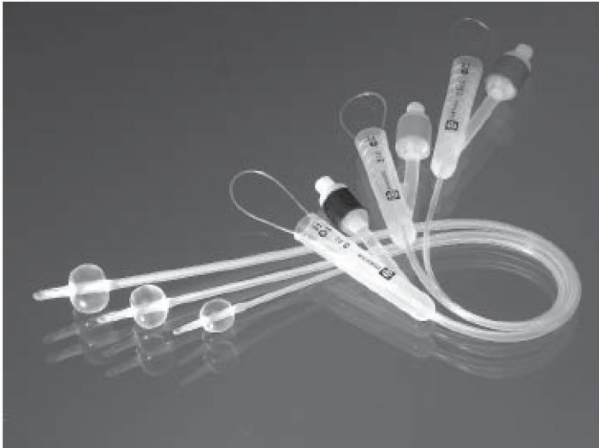

SPECIFICATIONS

2-Way Foley Balloon Catheter, Pediatric:

REF. No.	Size	Balloon	Length	Description
1821-0506	6 FR	1-1.5 cc/ml	330 mm	Stylet X-ray opaque line Catheter spigot
1821-0508	8 FR	2-3 cc/ml		
1821-0510	10 FR	3-5 cc/ml		
1821-0606	6 FR	1-1.5 cc/ml	330 mm	Metal Stylet X-ray opaque line Catheter spigot
1821-0608	8 FR	2-3 cc/ml		
1821-0610	10 FR	3-5 cc/ml		

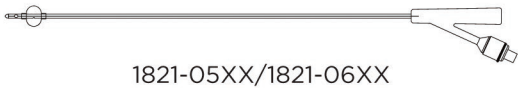

UROLOGICAL SERIES

SILICONE FOLEY BALLOON CATHETER


2-Way Foley Balloon Catheter, Standard, 5-10 cc/ml:

REF. No.	Size	Balloon	Length	Description
1822-0512	12 FR	5-10 cc/ml	420 mm	X-ray opaque line Catheter spigot
1822-0514	14 FR			
1822-0516	16 FR			
1822-0518	18 FR			
1822-0520	20 FR			
1822-0522	22 FR			
1822-0524	24 FR			
1822-0526	26 FR			




2-Way Foley Balloon Catheter, Standard, 20-30 cc/ml:

REF. No.	Size	Balloon	Length	Description
1823-0514	14 FR	20-30 cc/ml	420 mm	X-ray opaque line Catheter spigot
1823-0516	16 FR			
1823-0518	18 FR			
1823-0520	20 FR			
1823-0522	22 FR			
1823-0524	24 FR			
1823-0526	26 FR			

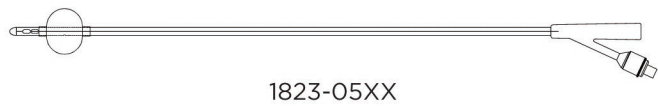

2-Way Foley Balloon Catheter, Female:

REF. No.	Size	Balloon	Length	Description
1824-0512	12 FR	5-10 cc/ml	220 mm	X-ray opaque line Catheter spigot
1824-0514	14 FR			
1824-0516	16 FR			
1824-0518	18 FR			
1824-0520	20 FR			
1824-0522	22 FR			
1824-0524	24 FR			
1824-0526	26 FR			

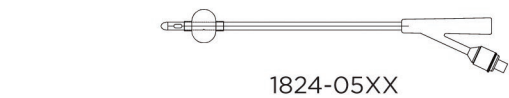


2-Way Foley Balloon Catheter, Tiemann Tip:

REF. No.	Size	Balloon	Length	Description
1856-0210	10FR	3-5 cc/ml	420 mm	Tiemann tip X-ray opaque line Catheter spigot
1856-0212	12 FR	5-10 cc/ml		
1856-0214	14 FR			
1856-0216	16 FR			
1856-0218	18 FR			
1856-0220	20 FR			
1856-0222	22 FR			
1856-0224	24 FR			



Catheter Spigot:

REF. No.	Size	Description
1820-0000	Universal	ABS Sterilized packaging 50 pcs/box

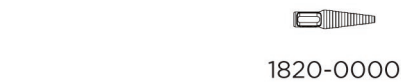

UROLOGICAL SERIES

SILICONE FOLEY BALLOON CATHETER

Fortune® Silicone 3-Way Foley Catheter Features:

- 100% medical grade silicone for superior biocompatibility.
- Transparent tube with X-ray opaque line.
- A soft and uniformly inflated balloon makes the tube sit well against the bladder.
- 3-Way design with an attached plug on the irrigation funnel.
- Tapered Central funnel could be connected with irrigation bags.
- Includes an individual sterilized packed Catheter Spigot.
- Sterilized double packaging.
- The tip of 3-way Foley Balloon Catheter with open hole, it could be used together with guide wire for difficulty in catheterization.


A Silicone 3-Way Foley Balloon



SPECIFICATIONS

3-Way Foley Balloon Catheter

REF. No.	Size	Balloon	Length	Description
1830-0514	14 FR	20-30 cc/ml	420 mm	Attached plug on flushing port X-ray opaque line Catheter spigot Open tip
1830-0516	16 FR			
1830-0518	18 FR			
1830-0520	20 FR			
1830-0522	22 FR			
1830-0524	24 FR			
1830-0526	26 FR			

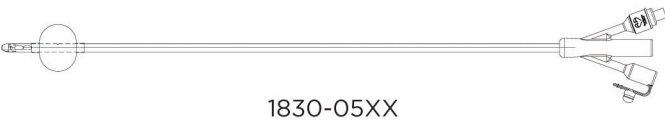

SILICONE FOLEY BALLOON CATHETER


Fortune® Foley Integrated Balloon Catheter:

Catheter with Integrated Balloon was designed to minimize the discomfort during catheterization procedure. Compare to the ordinary balloon, due to its thickness, creates a conjunction edge after adhered onto the tubing shaft. This conjunction edge between balloon and shaft will eventually increase patient discomfort during insertion and withdrawal.

Fortune® Integrated Balloon Catheter is invented to eliminate the uneven edge out of ordinary balloon. Patented balloon is integrated into the catheter shaft, creates “complete” smooth surface. No conjunction edge between balloon and shaft is formed, which minimizes the friction within the urethra and reduced patient discomfort.

Advantages:


- Patented balloon integrated into the shaft.
- Complete smooth surface.
- Reduce friction within the urethra
- Reduce trauma and secondary harm.
- Reduce discomfort during insertion and withdrawal.
- Reduce accumulation of urine sediment.

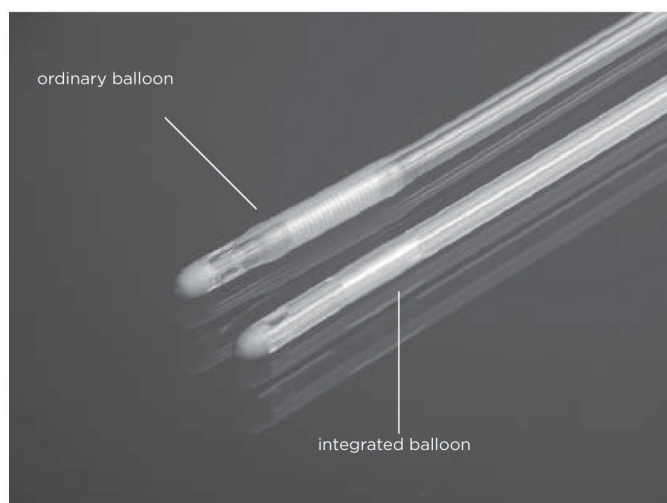

A Silicone Foley Integrated Balloon Catheter

B Silicone Foley Integrated Balloon Catheter, Tiemann

Tip

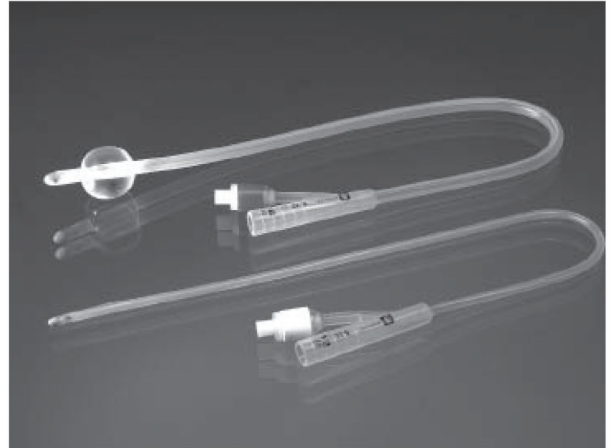

CE 0086


UROLOGICAL SERIES

SILICONE FOLEY BALLOON CATHETER

SPECIFICATIONS


Foley Integrated Balloon Catheter:

REF. No.	Size	Balloon	Length	Description
4822-0512	12 FR	5-10 cc/ml	420 mm	Integrated balloon X-ray opaque line Catheter spigot
4822-0514	14 FR			
4822-0516	16 FR			
4822-0518	18 FR			
4822-0520	20 FR			
4822-0522	22 FR			
4822-0524	24 FR			
4822-0526	26 FR			

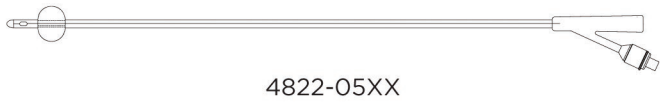

Foley Integrated Balloon Catheter, Tiemann Tip:


REF. No.	Size	Balloon	Length	Description
4856-0212	12 FR	5-10 cc/ml	420 mm	Tiemann tip Integrated balloon X-ray opaque line Catheter spigot
4856-0214	14 FR			
4856-0216	16 FR			
4856-0218	18 FR			
4856-0220	20 FR			
4856-0222	22 FR			
4856-0224	24 FR			

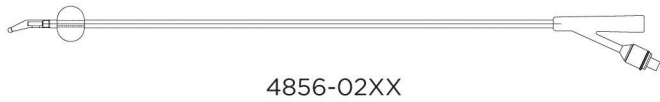

SILICONE FOLEY BALLOON CATHETER ,GROOVED



Fortune[®] Silicone 2 way Foley balloon grooved Catheter Features:

- 100% medical grade silicone for superior biocompatibility.
- Transparent medical grade silicone tube allows easy visual inspection and fluid observation.
- Grooved catheter was designed to minimize the discomfort during catheterization procedure. With ribs on outer profile of tubes, lubricity of catheter is increased, facilitating easy catheter removal.
- Soft and uniformly inflated balloon makes the tube sit well against the bladder.
- Includes an individual sterilized packed Catheter Spigot.
- Sterilized double packaging.

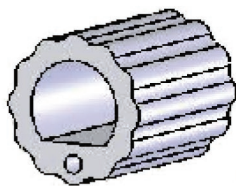

Grooved Tube

SPECIFICATIONS

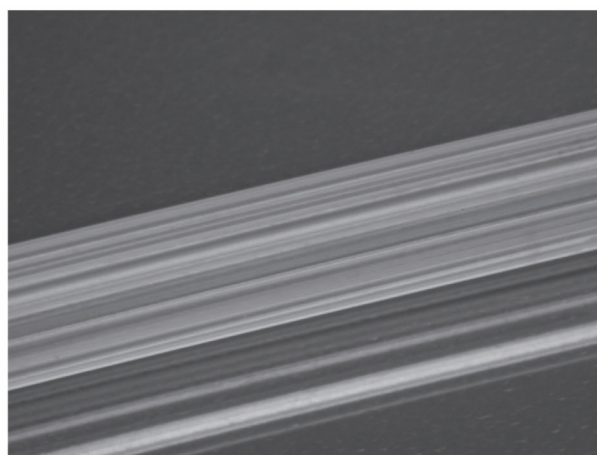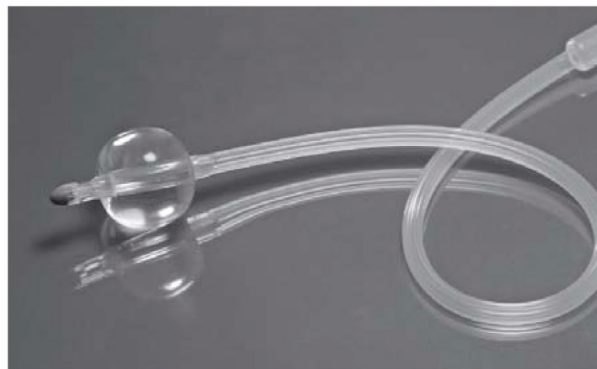

Silicone 2-way Foley balloon catheter, grooved tube

REF. No.	Size	Balloon	Length	Description
1827-0012	12 FR			
1827-0014	14 FR			
1827-0016	16 FR			
1827-0018	18 FR	5-10 cc/ml	420 mm	grooved tube, Catheter spigot
1827-0020	20 FR			
1827-0022	22 FR			
1827-0024	24 FR			



1827-0018



Tube Section

