

Veterinary ENT diagnostic instruments

RELIABLE, LONG-LASTING INSTRUMENTS

The HEINE veterinary range offers instruments of the highest quality and functionality as well as HEINE XHL Xenon halogen bulbs and innovative power sources.

HEINE VET instruments are subjected to ongoing quality assurance checks. They stand up to the demands of the busiest veterinary practice and assure optimum levels of diagnosis and treatment.



HEINE Veterinary ENT diagnostic instruments



Slit Illumination Head

BETA200 F.O.
VET OtoscopeBETA100
VET Otoscope

	Housing construction		
Metal	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
	Swivel lens – Magnification		
Glass/Magnification 2x	<input type="checkbox"/>		
Glass/Magnification 2,5x		<input type="checkbox"/>	<input type="checkbox"/>
	Illumination		
F.O.	<input type="checkbox"/>	<input type="checkbox"/>	
Direct			<input type="checkbox"/>
	Instrument size		
Professional	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
	Available power sources		
Battery 2.5V	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Rechargeable battery 3.5V	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
	Tips		
Long, closed specula	<input type="checkbox"/>		
Slit metal specula, open	<input type="checkbox"/>		
Short, closed specula	<input type="checkbox"/>		
UniSpec VET disposable tips	<input type="checkbox"/>		
Tip adaptor	<input type="checkbox"/>		
VET reusable specula		<input type="checkbox"/>	<input type="checkbox"/>
AllSpec disposable tips (with adaptor)			<input type="checkbox"/>
Insufflation port		<input type="checkbox"/>	<input type="checkbox"/>
Page	035	037	037

* with Adaptor

Slit Illumination Head and accessories



- ✦ **Solid metal construction.** Will last a lifetime.
- ✦ **Slit instrument head.** Simplifies the use of instruments.
- ✦ **Fiber optics illumination.** Bright, homogeneous light.
- ✦ **Multi-coated, scratch-resistant swivel lens with 1.6x magnification.** High-resolution image.
- ✦ **Simple open construction.** Easy to clean and disinfect with conventional solutions.
- ✦ **Double bayonet lock for specula.** Specula mounted firmly by double lock.

Slit Illumination Head	2.5V	3.5V
with swivel lens	G-001.21.301	G-002.21.301
XHL Xenon Halogen spare bulb	X-001.88.077	X-002.88.078
Swivel lens	G-000.21.209	

Long, closed specula

for Slit Illumination Head

For use with the Slit Illumination Head.



57 mm 65 mm 90 mm 150 mm
4 mm dia. 6 mm dia. 7 mm dia. 10 mm dia.

- ✦ **Metal construction, chrome-plated.** Easy to clean, disinfect and sterilize.
- ✦ **Double bayonet lock for specula.** Firm mounting of specula to head.

Closed specula for Slit Illumination Head		
57 mm long / 4 mm dia.	G-000.21.330	
65 mm long / 6 mm dia.	G-000.21.331	
90 mm long / 7 mm dia.	G-000.21.332	
150 mm long / 10 mm dia.	G-000.21.350	

Slit metal specula, open

for Slit Illumination Head

For use with the Slit Illumination Head.



65 mm long 90 mm long
6 mm dia. 7 mm dia.

- ✦ **High-quality chrome-plated metal construction.** Easy to clean, disinfect and sterilize.
- ✦ **Matt internal surfaces.** Reduced reflexes.
- ✦ **60° opening with rounded edges.** Ideal for instrumentation, atraumatic insertion.
- ✦ **Double bayonet fitting.** Positive fit.

Slit metal specula for Slit Illumination Head		
65 mm long / 6 mm dia.	G-000.21.341	
90 mm long / 7 mm dia.	G-000.21.342	

Short, closed specula

for Slit Illumination Head and BETA100 VET Otoscope



- ✦ **Made of high-density SANALON S with a smooth external surface.** Dust and germs are easily removed. Disinfection and autoclaving at up to 134 °C is possible.
- ✦ **Ergonomic shape.** Painless atraumatic insertion.
- ✦ **Integral metal ring with 2 bayonet locks.** Firm mounting of specula to head.
- ✦ **40 mm long.**

Short, closed SANALON S specula

for Slit Illumination Head and BETA100 VET Otoscope

1 set = 5 veterinary specula in SANALON S	G-000.21.316
40 mm long / 2.2 mm dia.	G-000.21.310
40 mm long / 2.8 mm dia.	G-000.21.311
40 mm long / 3.5 mm dia.	G-000.21.312
40 mm long / 4.5 mm dia.	G-000.21.313
40 mm long / 5.5 mm dia.	G-000.21.314

HEINE UniSpec® disposable ear specula

for BETA100 and with Adaptor for Slit Illumination Head



The adaptor illustrated [01] is necessary when using UniSpec specula with the Slit Illumination Head.

- ✦ **No cleaning necessary.** Saves time.
- ✦ **Hygienic and safe.** Eliminates cross-infection.
- ✦ **Packed in units of 1000.** Assures a competitive price.
- ✦ **Original HEINE design and manufacture.** Guaranteed good fit to the instrument head.

UniSpec disposable ear specula (grey)

Pack of 1000 UniSpec disposable specula 2.5 mm dia.	B-000.11.242
Pack of 1000 UniSpec disposable specula 4 mm dia.	B-000.11.241
Adaptor for use with Slit Illumination Head [01]	G-000.21.302

⊗: Single use only.

Spreadable speculum for Slit Illumination Head




This versatile instrument is widely used as a universal speculum for small animals.

- ✦ **For use with the Slit Illumination Head.**
- ✦ **Distal jaw opens up to 28 mm with 4 click-stops.** Useful for examination and therapy.
- ✦ **High-quality, chrome-plated metal construction.** Easy to clean, disinfect and sterilize.
- ✦ **Matt internal surfaces.** Avoids reflexes.
- ✦ **Oval branches with rounded edges.** Ideal for instrumentation, atraumatic.
- ✦ **2 bayonet locks.** Firm mounting of spreadable Speculum to head.

Spreadable speculum for Slit Illumination Head (95 mm long)	G-000.21.360
---	---------------------

HEINE BETA®200 Fiber Optic (F.O.) VET Otoscope

Veterinary Otoscope with Fiber Optic illumination

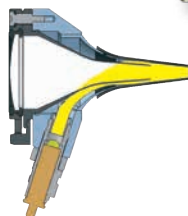

Fiber Optic (F.O.)
illumination

Solid, all-metal construction. A quality instrument for the Vet's office with Fiber Optic (F.O.) illumination. Insufflation port. The swivel window features an airtight seal.

- ✦ **Distal Fiber Optic illumination.** Perfectly even illumination of auditory canal and tympanum.
- ✦ **Large, scratch-resistant glass viewing window with 2.5x magnification.** Unobstructed view, distortion-free image.
- ✦ **Swivelling viewing window built into instrument.** Useful for instrumentation, cannot be mislaid.
- ✦ **Insufflation port.** Reliable test of tympanic mobility.
- ✦ **Durable metal body, matt chrome-plated.** Unbreakable, long-life.
- ✦ **XHL Xenon Halogen Technology.** Extremely bright, white light.
- ✦ **For use with reusable SANALON S VET Specula (G-000.21.213).**
- ✦ **Instrument head matt-black inside.** Eliminates reflexes.

BETA200 VET Otoscope	2.5V	3.5V
without specula, without handle	G-001.21.250	G-002.21.250
XHL Xenon Halogen spare bulb	X-001.88.077	X-002.88.078

HEINE BETA®100 VET Otoscope

Veterinary Otoscope with direct illumination



Direct Xenon Halogen
illumination

Otoscope with all-metal construction and direct Xenon Halogen illumination. Suitable for use with short, closed specula. With the use of the VET adaptor, it is possible to fit reusable VET specula for ear examinations.

- ✦ **All-metal housing.** Exceptionally-tough and long-lasting.
- ✦ **XHL Xenon Halogen Technology.** Extremely bright, white light.
- ✦ **Large, scratch-resistant glass viewing window with 2.5x magnification.** Unobstructed view, distortion-free image.
- ✦ **Swivelling viewing window built into instrument.** Useful for instrumentation, cannot be mislaid.
- ✦ **Insufflation port.** Reliable test of tympanic mobility.
- ✦ **Can be used with short, closed specula G-000.21.316 (no adaptor required).**
- ✦ **Versatile, useful for very small animals.**
- ✦ **Instrument head matt-black inside.** Eliminates reflexes.

BETA100 VET Otoscope	2.5V	3.5V
without specula, without handle	G-001.21.260	G-002.21.260
XHL Xenon Halogen spare bulb	X-001.88.037	X-002.88.049
Adaptor for reusable SANALON S VET specula [01]	G-000.21.214	
Adaptor for AllSpec tips [02]	B-000.11.306	

[01]

[02]

Reusable VET specula

for BETA200 or BETA100 with adaptor



- ✶ **Made of high-density SANALON S with a smooth external surface.** Dust and germs are easily removed. Disinfection and autoclaving at up to 134 °C is possible.
- ✶ **Ergonomic shape.** Atraumatic.
- ✶ **Convenient fixation system.** Positive fit to otoscope head.

Reusable SANALON S specula

1 set = 3 reusable SANALON S specula	G-000.21.213
61 mm long / 4 mm dia.	G-000.21.210
61 mm long / 6 mm dia.	G-000.21.211
61 mm long / 9 mm dia.	G-000.21.212

HEINE AllSpec® Disposable tips


for BETA200 or BETA100 with adaptor



- ✶ **For single use.** Hygienic, no cross-infection.
- ✶ **Original HEINE design.** Firm fit to the otoscope head.
- ✶ **Assured quality.** No sharp edges.

AllSpec Disposable tips (grey)

Pack of 1000 pcs. 2.5 mm dia. [03]	B-000.11.128
Pack of 1000 pcs. 4 mm dia. [03]	B-000.11.127
Case of 10 dispenser packs of 250 tips, 2.5 mm dia [01]	B-000.11.151
Case of 10 dispenser packs of 250 tips, 4 mm dia [01]	B-000.11.150
AllSpec Tip-dispenser [02]	B-000.11.148

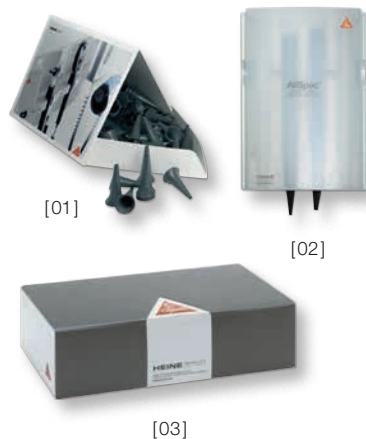

AllSpec Tip-dispenser [02]

The right size, visible and instantly available. A clean new tip for each and every patient. Contents: 50 pcs. Ø 2.5 mm and 50 pcs. Ø 4 mm.

⊗: Single use only.

Insufflation Bulb

For BETA200VET and BETA100VET Otoscope.



- ✶ **For pneumatic testing.**

Insufflation bulb	B-000.11.240
-------------------	---------------------

HEINE Veterinary Diagnostic Sets



Set complete with: **Slit illumination head**

1 set long, closed specula (57x4 mm, 65x6 mm, 90x7 mm)

1 spare bulb

hard case



Veterinary Diagnostic Sets	2.5V	3.5V (NiMH)	3.5V 
BETA battery handle	G-112.10.118		
BETA R rechargeable handle for mains socket		G-112.20.376	G-112.27.376
BETA TR rechargeable handle with plug-in transformer		G-112.20.384	G-112.27.384
BETA NT rechargeable handle and NT300 charger		G-112.20.420	
 BETA L rechargeable handle and NT300 charger			G-112.29.420



Set complete with: **Slit illumination head, BETA200 Ophthalmoscope**

1 set long, closed specula (57x4 mm, 65x6 mm, 90x7 mm)

1 spare bulb for each instrument

hard case





Veterinary Diagnostic Sets	2.5V	3.5V (NiMH)	3.5V 
BETA battery handle	G-148.10.118		
BETA R rechargeable handle for mains socket		G-148.20.376	G-148.27.376
BETA TR rechargeable handle with plug-in transformer		G-148.20.384	G-148.27.384
BETA NT rechargeable handle and NT300 charger		G-148.20.420	
 BETA L rechargeable handle and NT300 charger			G-148.29.420

: Rechargeable handle with charge status display in the bottom insert.

HEINE Veterinary Diagnostic Sets





Set complete with: **Slit illumination head, BETA200 Ophthalmoscope**
 5 short, closed SANALON S specula
 10 UniSpec disposable ear specula with adaptor
 1 set = 4 long, closed specula (57 x 4 mm, 65 x 6 mm, 90 x 7 mm, 150 x 10 mm)
 2 slit specula
 spreadable speculum
 1 spare bulb for each instrument
 hard case

Veterinary Diagnostic Sets	2.5V	3.5V (NiMH)	3.5V 
BETA battery handle	G-150.10.118		
BETA R rechargeable handle for mains socket		G-150.20.376	G-150.27.376
BETA TR rechargeable handle with plug-in transformer		G-150.20.384	G-150.27.384
BETA NT rechargeable handle and NT300 charger		G-150.20.420	
 BETA L rechargeable handle and NT300 charger			G-150.29.420



Set complete with: **BETA200 Ophthalmoscope, BETA200 VET Otoscope**
 3 reusable SANALON S VET specula
 1 spare bulb for each instrument
 zipper case

Veterinary Diagnostic Sets	2.5V	3.5V (NiMH)	3.5V 
BETA battery handle	G-151.10.118		
BETA R rechargeable handle for mains socket		G-151.20.376	G-151.27.376
BETA TR rechargeable handle with plug-in transformer		G-151.20.384	G-151.27.384
BETA NT rechargeable handle and NT300 charger		G-151.20.420	
 BETA L rechargeable handle and NT300 charger			G-151.29.420

: Rechargeable handle with charge status display in the bottom insert.