

Heat for infusions
and blood products



The Barkey *S-line* smart – simple – safe

Advanced technology: One-touch-control

- The Barkey s-line warms blood products and infusions up to the patient.
- The electrically heated and single-sided opening heating device is capable of accepting standard IV-sets.
- The Barkey s-line is simple to operate in all medical departments.
- It is extremely easy to clean due to its smooth surface.
- Its large and flexible application areas at low to medium flow rates result in an optimum alternative to cost-intensive blood heaters with single-use articles.


Avoidance of hypothermia

Warmth is the key to your patients' well-being! With Barkey's inexpensive s-line devices you can apply warmth which eases treatment for all patients and accelerates the healing process. The result: care personnel has less work and treatment costs are reduced.

Made-to-order technology for medicine

Summary of further benefits

- Integrated temperature sensors
- guarantee best warming and highest safety.
- Complete casing of IV-sets and therefore optimum warming
- No additional consumer articles
- Warming of blood return feeds or of dialysis fluids in combination with devices for replacement kidney therapy
- Warming of enteral and/or parenteral nutrition
- Operational safety thanks to preset target temperature value
- Display of actual temperature of the heating sleeve
- Available in a range of diameters
- Smooth surface, therefore simple to wash and disinfect
- Use of exclusively bio-compatible materials
- Optical and acoustic alarm signal in fault situation

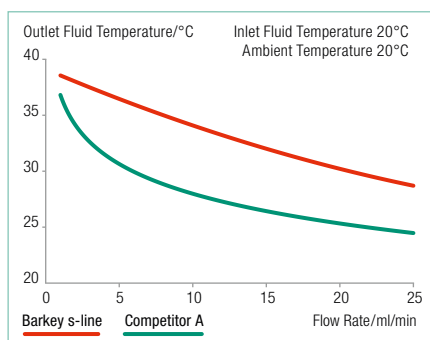

Exclusive innovative technological solutions at the highest quality standards and specialised benefits in practice - this is the basis for the success of the holistic Barkey heating concept. Barkey has earned its excellent reputation with exclusive engineering know-how in medical technology over the last 25 years. The proverbial quality and practical suit-

ability of the devices and systems "made in Germany" has established itself in a wide range of applications.

Real innovations and problem solutions which provide the user with specific benefits can only be developed by people who keep their ears and eyes open, who talk to specialists in practice and recognise their requirements in good time.

For this reason we would like to invite you to carry out an active dialogue with us: Have you any questions about devices and systems in the Barkey heating concept?

Which specialist questions can we answer for you? What suggestions have you got for us?

Tel.: +49 (5202) 9801-0

Fax: +49 (5202) 9801-99

e-mail: info@barkey.de

Type	Barkey s-line
Article number	213.10285
Temperature setting	39°C setpoint temperature ex-works
Heating profile length	L: 1.500 mm
Internal diameters available	4.55 mm 6.55 mm
Suitable for infusion piping diameters	of 4.0 to 7.0 mm
Weight	approx. 750 g
Mains supply	100-240V, 50/60Hz
Connection	24 volt DC
Power consumption	60 Watt
Connecting cable	3 m
Safety class	I protective insulation Degree of protection BF external and internal application to patients with insulated application parts
Humidity protection	IP X2
Certification	CE 0123
Barkey is certified to	ISO 9001 ISO 13485
Made in Germany	



The accessories shown are not included in delivery.

Sold by:

www.barkey.de