



**Cricath<sup>®</sup>, the smart choice  
for cricothyrotomy**






### **Cricath<sup>®</sup>, the smart choice:**

- › small lumen – 2 mm ID
- › of flexible, high grade Teflon<sup>®</sup>
- › and an innovative tilting mechanism
- › to reduce the risk of kinking

Cricath<sup>®</sup> is a flexible needle cricothyrotomy catheter that enables a quick access to the trachea. This application asks for the best materials. The lumen of Cricath<sup>®</sup> is made of a special, high grade type of Teflon<sup>®</sup> (PFA). The patented tilting mechanism enables a proper fit of the lumen into the trachea, also in case of an abnormal neck anatomy. Together with the elevated fixation point of the lumen and the enlarged, smoothly curved opening in the flange around the lumen, this provides for additional flexibility, reducing the risk of kinking. With the standard catheter over needle procedure for cricothyrotomy, Cricath<sup>®</sup> can be placed by insertion through the cricoid membrane.

cricath<sup>®</sup>, the smart choice

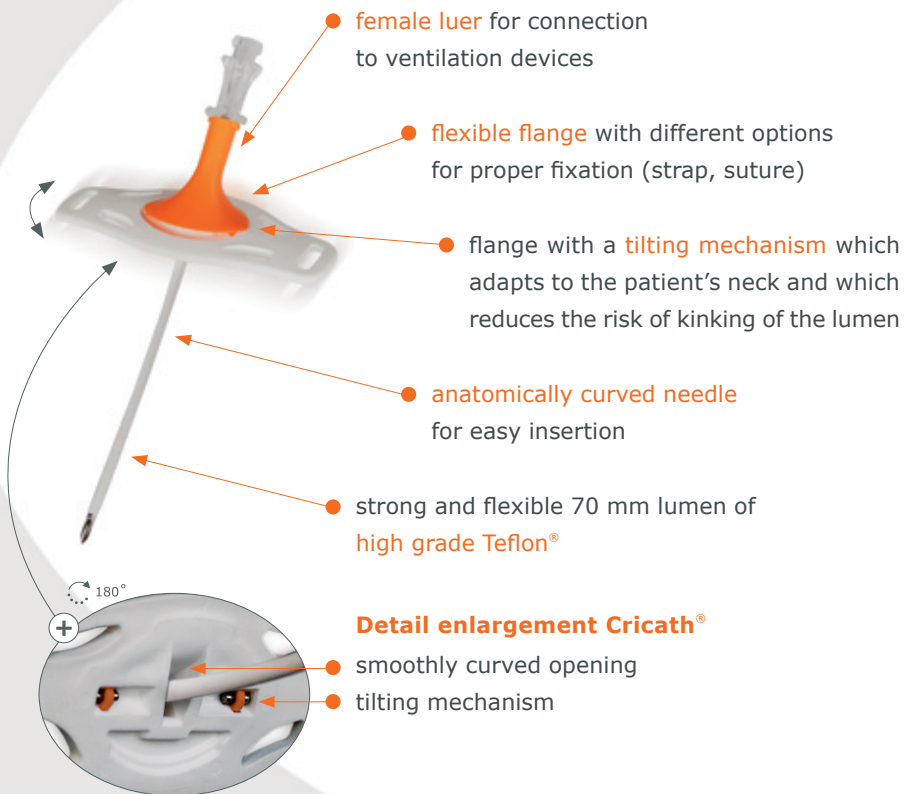

*Cricath® is protected by a patent application and a design registration.*



**quick  
access to  
the trachea**

ce for cricothyrotomy



*Paola Vasters  
product manager*

*Dirk van Asseldonk  
CEO*

*Caroline Abeln  
customer services*



### **Development**

Cricath<sup>®</sup> is developed in the Netherlands by Dolphys Medical, the same company that developed Ventrain<sup>®</sup>. For more information about Cricath<sup>®</sup> and our other products please visit our website or contact us directly.

### **Ordering information Cricath<sup>®</sup>**

- packed with syringe and neck strap
- single use
- sterile
- art. code 2618400

### **Website**

[www.dolphys.com/cricath](http://www.dolphys.com/cricath)

#### **Your distributor**

#### **Dolphys Medical B.V.**

**Info** Paola Vasters  
Caroline Abeln

**P** De Lismortel 31  
5612 AR Eindhoven  
The Netherlands

**T** +31 (0)40 751 60 20  
**E** [cricath@dolphys.com](mailto:cricath@dolphys.com)  
**I** [www.dolphys.com](http://www.dolphys.com)