


## TRACOE *vario*

January 2012

TRACOE *vario extract* with an integrated suction line for subglottic suctioning is now also available in extra-long as REF 471!

The TRACOE *vario extract XL* with adjustable neck flange is designed to minimise the risk of airway infections in obese patients.

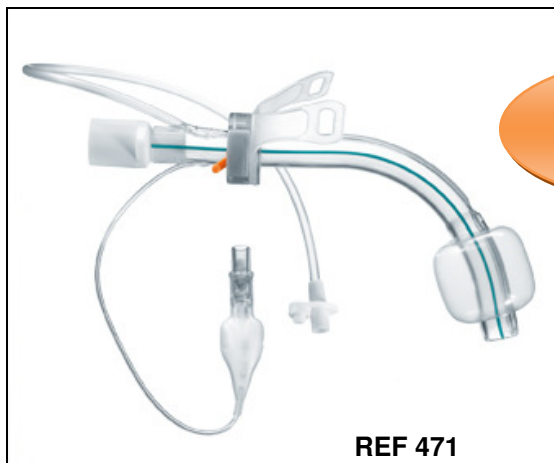

**NEW!**

### Function:

The tracheostomy tube has a suction line which is integrated into the tube wall with an opening right above the cuff. Through this opening of the suction line the accumulated secretion can be removed.

### Advantages:

- Extra-long tracheostomy tube with subglottic suctioning for obese patients.
- The risk of transferring bacterially contaminated secretions into the lower respiratory tract is reduced and therefore also the risk of developing aspiration pneumonia (e.g. VAP).
- Possible therapy support for patients with dysphagia.
- Through the integration of the suction line into the tube wall, pressure marks in the tracheostoma will be avoided.
- Through the variably adjustable neck flange with rotating wings a universal adjustment to any patient's neck anatomy is possible (e.g. thick neck, postoperative oedema, asymmetrical necks).

### Dimensions REF 471

Size	ID (mm)	OD (mm)	Length TL* (mm)	Length AB** (mm)
07	7.0	9.7	104	110
08	8.0	11.2	123	128
09	9.0	12.3	126	134
10	10.0	13.7	135	146

\*TL: length along centre line - \*\* AB: length along outer bend

For further information please visit our homepage: [www.tracoe.com](http://www.tracoe.com)